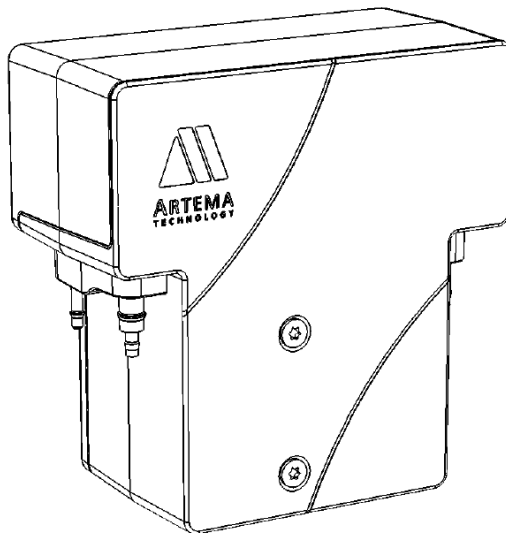


# AION Rhodium® CO<sub>2</sub> Analyzer

## Product Specification



*Interfaces to optional*



### INTRODUCTION

AION Rhodium® is a highly integrated CO<sub>2</sub> analyzer designed specifically for the OEM market. Even though it is not classified as a finished medical device it fulfils all relevant regulatory requirements. This, combined with a high degree of integrated parameters provides for an uncomplicated integration into host systems. The AION Rhodium® interfaces to optional oxygen sensors, paramagnetic or galvanic and can also be used together with the SPIRIT™ Respiratory Mechanics Analyzer for a full set of respiratory parameters.

**INTENDED USE**

AION Rhodium® is a highly integrated gas analyzer intended to continuously analyze CO<sub>2</sub> concentration and breath rate during intensive care and anesthesia. When used with an optional oxygen sensor, paramagnetic or galvanic, the analyzer also provides O<sub>2</sub> concentration data. All parameters are transmitted to a host instrument over a digital data interface. It is a compact, low-power and high-performance sub-unit intended to be integrated in other manufacturer's medical devices.

*Indications:*

- The respiratory gas measurement information provided by the analyzer is intended to be used by trained and authorized healthcare professionals only.
- The analyzer is only intended to be integrated into host devices by engineering professionals with relevant knowledge of medical device development and manufacturing. The host device manufacturer is responsible for the correct integration of the analyzer and handling of any risks that may arise.
- For accurate and safe function, the AION Rhodium® shall be used only together with Artema Technology® subsystems, accessories and disposables.
- The host instrument must fulfill applicable national requirements for Medical Devices.

*Contraindications:*

- The analyzer is not intended for use with flammable anesthetic agents.
- The analyzer is not intended to be used in outdoor transport applications such as cars or aircrafts.
- The analyzer is not intended to be used as the only means of monitoring a patient and it shall always be used in combination with other vital sign monitoring devices and/or professional human judgments of patient condition.
- **WARNING:** No modification of the analyzer is allowed.

NORMAL RANGE GENERAL SPECIFICATIONS			
<b>Technology</b>	NDIR type, diverting (side-stream) gas analyzer with optional external paramagnetic or galvanic oxygen sensor. Automatic ambient pressure and temperature compensation. Host supplied O <sub>2</sub> <sup>1)</sup> , N <sub>2</sub> O and agent compensation input.		
<b>Oxygen Sensor</b>	<ul style="list-style-type: none"> <li>• Optional external Hummingbird Paracube® Sprint Paramagnetic Oxygen sensor</li> <li>• Optional external OXIMA™ Galvanic Oxygen Sensor</li> </ul>		
<b>Operating modes</b>	<ul style="list-style-type: none"> <li>• Startup</li> <li>• ISO accuracy</li> <li>• Full accuracy</li> </ul>		
<b>Measured gases</b>	CO <sub>2</sub> , and O <sub>2</sub> <sup>2)</sup>		
<b>Measured parameters</b>	<ul style="list-style-type: none"> <li>• Momentary gas concentration</li> <li>• Inspired and expired concentrations of all measured gases</li> <li>• Breath rate</li> </ul>		
<b>Resolution</b>	CO <sub>2</sub> : 0,01%; O <sub>2</sub> <sup>2)</sup> : 0,1%		
<b>CO<sub>2</sub> Zero threshold</b>	Gas concentration set to zero if level below: <ul style="list-style-type: none"> <li>• in full accuracy mode: 0.1 %</li> <li>• in ISO accuracy mode: 0.3%</li> </ul> for more than 3 s		
<b>Initialization time</b>	Performs a self-test and enters standby mode within 10s where no data frames will be transmitted after hardware reset or a Reset command.		
<b>System rise times (t<sub>10-90</sub> %)</b> • 70 ml/min using DRYLINE™ II • 120 ml/min using DRYLINE™ II	<b>Gas component</b>	<b>@ 70 ml/min</b>	<b>@ 120 ml/min</b>
	CO <sub>2</sub>	200 ms	200 ms
	O <sub>2</sub> Paramagnetic <sup>2)</sup> 15 to 21% (21 to 60%)	/	600 ms (800 ms)
	O <sub>2</sub> OXIMA™ <sup>2)</sup> 15 to 21% (21 to 60%)	/	500 ms (600 ms)
<b>System<sup>3)</sup> delay time (t<sub>0-10</sub> %)</b>	<4 seconds		
<b>Patient category</b>	Adult, pediatric and neonatal		

**Notes:**

1. Automatic O<sub>2</sub> compensation if used with optional oxygen sensor
2. With optional oxygen sensor.
3. For a complete system with Artema Technology gas sampling accessories.

CO <sub>2</sub> MEASUREMENT SPECIFICATIONS <sup>1)</sup>	
Full accuracy	±0,263 % <sub>abs</sub> (±2 mmHg)
ISO accuracy <sup>2)</sup>	±(0,43 % <sub>abs</sub> + 8 % <sub>rel</sub> )

**Notes:**

1. Inaccuracy specifications include stability and drift.
2. Corresponding to requirements in ISO 80601-2-55:2011

CO <sub>2</sub> ACCURACY MODE APPLICATION			
The AION® Rhodium provides full accuracy specification (white areas in the graphs below) in a majority of clinical cases and environmental conditions. During more extreme conditions, clinical or environmental, the AION Rhodium® provides ISO accuracy (light grey areas below).			
Gas range	<div>CO<sub>2</sub> level 0 – 6 %: Full</div> <div>CO<sub>2</sub> level 6 – 20 %: ISO</div>		
Time from startup	<div>&lt;10 s Startup</div> <div>10 s – 5 min: ISO</div> <div>&gt; 5 min: Full</div>		
Ambient temperature	<div>10 – 15 °C : ISO</div> <div>15 – 30 °C: Full</div> <div>30 – 55 °C: ISO</div>		
Ambient pressure	<div>660 – 800 hPa: ISO</div> <div>800 – 1100 hPa: Full</div> <div>1100 – 1200 hPa ISO</div>		
End tidal value accuracy vs. respiration rate <i>For I:E ratio 1:1</i>	<div>0 – 60 bpm: Full<sup>1)</sup></div> <div>60 – 150 bpm: ISO</div>		

**Notes:**

1. For CO<sub>2</sub> range 0 to 6%

O <sub>2</sub> MEASUREMENT SPECIFICATIONS <sup>1), 2)</sup>				
Gas	Gas level [%]	Inaccuracy [% <sub>ABS</sub> ]	Interference [% <sub>ABS</sub> ]	Resp. rate limit for accurately resolved EtO <sub>2</sub> <sup>3)</sup>
O <sub>2</sub> Paracube® Sprint	0 – 25 25 – 80 80 – 100	±1 ±2 ±3	CO <sub>2</sub> 0.2 N <sub>2</sub> O 0.2 Any agent 1.0	30 bpm
O <sub>2</sub> OXIMA™	0 – 40 40 – 60 60 – 80 80 – 100	±(1%abs +1%rel) ±(1%abs +2%rel) ±(1%abs +3%rel) ±(1%abs +4%rel)	<0.3	60 bpm

**Notes:**

1. With optional oxygen sensor
2. Inaccuracy specifications include stability and drift.
3. At 120 ml/min using Artema Technology® neonate gas sampling accessories

GAS CONTAMINANTS INTERFERENCE			
Contaminant	Interference [% <sub>ABS</sub> ]		
	CO <sub>2</sub>	O <sub>2</sub> Paracube® Sprint (optional)	O <sub>2</sub> OXIMA™ (optional)
< 100 % Xenon	0.1	0.5	0.3
< 50 % He	0.1	0.5	0.3
Metered dose inhaler propellants	Unspec.	0.5	Unspecified
< 0.1 % Ethanol	0	0.5	0.3
Saturated Isopropanol vapour	0.1	0.5	Unspecified
< 1 % Acetone	0.1	0.5	0.3
< 1 % Methane	0.1	0.5	0.3

PNEUMATIC SPECIFICATIONS	
<b>Technology</b>	Side-stream gas sampling
<b>Pneumatic modes</b>	<ul style="list-style-type: none"> <li>Room air reference measurement: Automatic</li> <li>Sampling system purge: Automatic</li> </ul>
<b>Pump</b>	Flow controlled membrane type
<b>Gas flow</b>	Gas sample and purge flows
<b>Gas sampling accessories:</b>	DRYLINE™ range of airway adapters, gas sampling lines and water traps
<b>Gas sampling rate<sup>1)</sup></b>	CO <sub>2</sub> only mode: 70 ml/min CO <sub>2</sub> + O <sub>2</sub> mode: 120 ml/min Flow accuracy ±10 ml/min or 10%, whichever is greatest.
<b>Gas system leakage</b>	<0.1% of set sampling flow
<b>Occlusion alarm</b>	Actual flow <40 ml/min
<b>Reference measurement interval in full accuracy mode<sup>2)</sup></b>	Automatic when gas measurement bench temperature change is > 1°C or time since last ref. measurement is > 4 hours.
<b>Reference measurement gas requirement</b>	Air (room or other) < 800 ppm CO <sub>2</sub>
<b>Reference measurement duration</b>	Typical: 4 s Maximum: 9 s <sup>3)</sup> Synchronized with inspiratory section of breathing cycle <sup>4)</sup> .
<b>Room air gas valve check</b>	Automatic, once every 24 h
<b>Purge cycle</b>	Automatic when occlusion detected
<b>Change water trap alarm</b>	Actual flow <75% of set flow
<b>Pressure difference<sup>5)</sup></b>	Pressure highest at gas sampling point: $P_{\text{SAMPLING POINT}} - P_{\text{EVAC}} : <30 \text{ hPa}$ Pressure highest at gas return point: $P_{\text{EVAC}} - P_{\text{SAMPLING POINT}}$ (normal operation): <40 hPa $P_{\text{EVAC}} - P_{\text{SAMPLING POINT}}$ (reference measurement): <120 hPa

**Notes:**

1. AION™ measures volumetric flow at actual barometric pressure, normalized to room air at 21 °C and 0% RH. The use of gas mixtures other than room air for flow calibration may cause flow measurement errors.
2. Reference measurement interval in ISO Accuracy Mode > 30 seconds.
3. Longer at high O<sub>2</sub> concentrations when using OXIMA™ galvanic oxygen sensor.
4. Synchronized with expiratory plateau of capnogram when reference valve function test performed.
5. For a complete system with Artema Technology gas sampling accessories.

GAS DATA OUTPUT	
<b>Breath detection</b>	Adaptive threshold (>1% <sub>ABS</sub> change in CO <sub>2</sub> concentration)
<b>Respiration rate</b>	2 – 60 bpm, accuracy ±1 bpm 60 – 150 bpm, accuracy ±2 bpm
<b>Fi and ET values</b>	CO <sub>2</sub> and O <sub>2</sub> <sup>1)</sup>
<b>Waveforms</b>	Up to five simultaneously
<b>Units for Measure<sup>2)</sup></b>	Gas data can be reported in [%], [hPa] or [Torr], separately selectable for CO <sub>2</sub> and O <sub>2</sub> <sup>1)</sup> .
<b>Status</b>	Atmospheric pressure, Sampling flow, Sensor head temperature and an extensive set of diagnostic parameters.
<b>Flags</b>	Apnea, No water trap, Water trap type, Change Water trap, Occlusion, Purge, Hardware errors.

**Notes:**

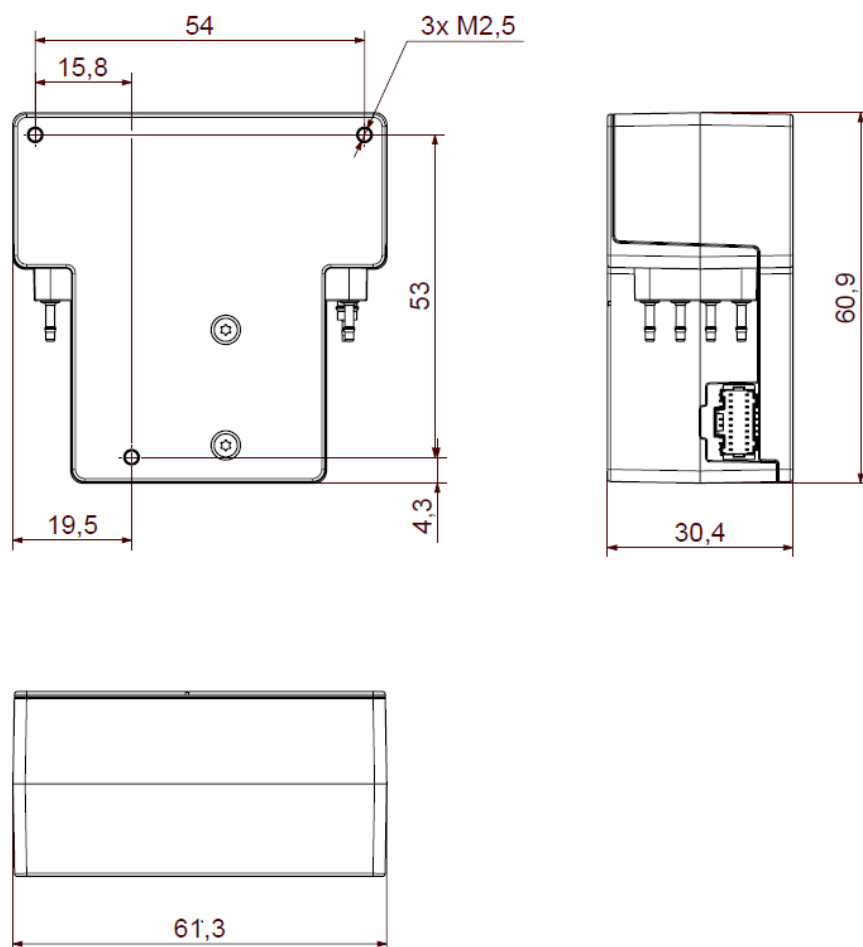
- 1 With optional oxygen sensor.
- 2 All gas readings are normally referenced to dry gas conditions, ambient room temperature and atmospheric pressure (ATPD). The gas readings may instead be referenced to saturated breathing gas at body temperature (BTSP) by sending a command to AION™ via the communication interface.


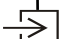
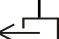
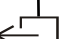
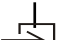

ELECTRICAL SPECIFICATIONS	
<b>Supply voltage</b>	11 – 30 V <sub>DC</sub> (max ripple 500 mV <sub>p-p</sub> )
<b>Power consumption<sup>1)</sup></b>	1.0 W Typical 1.8 W Maximum
<b>Electrical interface</b>	One polarized 15-pin connector for all electrical interfaces.
<b>Communications interface</b>	AION™ standard communications protocol; RS-232, 10ms data interval, 115.2 kBaud.
<b>Data sample rate</b>	50Hz. Data presentation is 100 Hz, every second data point is interpolated.

**Notes:**

- 1 When powered at 12 VDC. Add maximum 0.3 W when powered at higher voltage

## MECHANICAL SPECIFICATIONS



<b>Weight</b>	265 g
<b>Color</b>	Silver
<b>Mechanical interface</b>	3 x threaded holes on one side for M2.5 bolts.
<b>Pneumatic interface</b>	<p>6 x Stainless steel fittings:</p> <ul style="list-style-type: none"> <li> Secondary flow inlet from water trap</li> <li> Gas sample inlet from water trap</li> <li> Evacuation outlet</li> <li> External oxygen sensor outlet</li> <li> External oxygen sensor return inlet</li> <li> External reference gas inlet</li> </ul>



<b>ENVIRONMENTAL SPECIFICATIONS</b>	
The device is not intended to be used outside of the specified environmental limits	
<b>Storage temperature</b>	-40 to 70 °C.
<b>Storage humidity<sup>1)</sup></b>	5 – 100% RH, condensing
<b>Storage pressure</b>	50 – 120 kPa
<b>Operating temperature<sup>2)</sup></b>	10 – 55 °C. Full accuracy within 15 – 30 °C
<b>Operating ambient humidity</b>	10 – 95% RH, non-condensing
<b>Operating ambient pressure</b>	660 – 1200 hPa. Full accuracy within 800 – 1100 hPa
<b>Operating ambient CO<sub>2</sub> concentration</b>	< 800 ppm
<b>Maximum positive pressure in patient breathing circuit</b>	120 hPa (during lung recruitment)
<b>Maximum negative pressure in patient breathing circuit</b>	-310 hPa (during closed suctioning)
<b>EMC compatibility</b>	Supports EMC requirements for life-supporting ME equipment as per IEC 60601-1-2:2007. Final EMC compliance is the responsibility of the host system integrator.
<b>MR compatibility</b>	The device is not intended for use in a MR environment without extra protective measures.

**Notes:**

- 1 After storage in a condensing atmosphere, the unit shall before use be kept for more than 24 h in an environment equivalent to the operating atmosphere.
- 2 Refers to the temperature of the cover.








REGULATORY COMPLIANCE (APPLICABLE PARTS)	
<b>ISO 80601-2-55:2011</b>	Medical electrical equipment – Particular requirements for the basic safety and essential performance of respiratory gas monitors.
<b>IEC 60601-1:2005</b>	Medical electrical equipment – General requirements for basic safety and essential performance.
<b>IEC 60601-1-2:2007</b>	Medical Electrical Equipment – Part 1-2: General Requirements for Safety – Collateral Standard Electromagnetic compatibility – Requirements and tests.
<b>IEC 62304:2006</b>	Medical device software – Software life cycle processes.
<b>RoHS 2011/65/EU</b>	Directive on the restriction of the use of certain hazardous substances in electrical and electronic equipment.

**Notes:**

The above standards refer to several sub-standards not listed in this document. It is the responsibility of the customer to ensure that the host instrument meets the requirements in the international standards and regulatory requirements relevant for the market where the host instrument is sold.

REGULATORY DESIGNATIONS	
<b>Essential performance</b>	<i>"Measurement accuracy" only, "Gas reading alarm conditions" and "Technical alarm conditions" are handled by the host device.</i>
<b>Protection against electric shock</b>	Non-electrical connection to patient. Protection against electrical shock is handled by host device.
<b>Oxygen rich environments</b>	The device is suitable for use in an oxygen rich environment.
<b>Mode of operation</b>	Continuous
<b>Ingress of water and particulate matter</b>	As the AION Rhodium® is a component and no classification according to IEC 60529 is made.
<b>Protection against Hazards of Explosion</b>	Not protected (ordinary)

PACKAGE INFORMATION	
<b>Packaging</b>	<ul style="list-style-type: none"> <li>• Transparent Static Shielding Bag</li> <li>• Die cut dissipative foam fittings</li> <li>• White mini-well box</li> </ul>
<b>Package size</b>	118 × 81 × 32 [mm]
<b>Package weight</b>	< 360 g including product.
<b>Package durability</b>	Meets "basic values" in ISO 4180: 2009: Storage, Road, Rail, Water and Air transport
<b>Package labeling</b>	Includes Product name, Part number, Firmware version, Serial number, Artema Technology logo, Mindray Medical Sweden AB name.

PACKAGE MARKING	
	Catalogue number
	Serial number
	Date of manufacture
	Manufacturer
	Storage humidity limitation
	Storage atmospheric pressure limitation
	Storage temperature limitation

ORDERING INFORMATION	
Designation	Part no.
AION Rhodium® CO <sub>2</sub> Analyzer <i>For use with external galvanic, or no oxygen sensor</i>	115-027548-00
AION Rhodium® CO <sub>2</sub> + O <sub>2</sub> Analyzer <i>AION Rhodium® and Hummingbird Paracube® SPRINT preassembled and calibrated.</i>	115-030104-00

**mindray**  
SWEDEN

**Mindray Medical Sweden AB**

Esplanaden 3E · SE – 172 67 Sundbyberg · Sweden

Tel +46 (0)8 555 54 100

[www.mindray.com/artema](http://www.mindray.com/artema) · [artemainfo@mindray.com](mailto:artemainfo@mindray.com)