

## Model FS6122

SIARGO MEMS FLOW SENSING PRODUCTS

MEMS Sensor Fusion for CPAP Applications

VA.5

















# MEMS Sensor Fusion for CPAP Applications

FS6122 Series

## **User Manual**

Document No. 10-2016-S8 EN

ssue date: 2016.10 Revision: VA.5

## Siargo Ltd.

3100 De La Cruz Boulevard, Suite 210, Santa Clara, CA 95054 USA

Tel: +1(408)969.0368 Email: info@siargo.com

## © Copyright 2016 by Siargo Ltd.

Siargo Ltd. and its subsidiaries reserve the rights to change the specifications and/or descriptions without prior notice. For further information and updates, please visit: www.Siargo.com

## **Medical Respiratory Sensor**

#### Model FS6122

# Siargo Ltd.

#### 1. Features

- Specially designed for medical CPAP applications
- Direct measurement of mass flow rate and gauge pressure
- Small dead volume with compact form factor
- Digital and analog linear output with fast response time
- Low power consumption



#### 2. Sensor Performance

All data unless otherwise noted apply for calibration conditions: air,  $20\,^{\circ}$ C,  $101.325\,kPa$  absolute pressure, horizontal mounting.

Spec	Vaule	Unit	
Flow Sensor			
Flow Range <sup>1</sup>	-250 ~ <b>+</b> 250	SLPM	
Accuracy (Total Error Band)	± (2.5 + 0.5FS)	%	
Output <sup>2</sup>	Linear, Analog/l²C		
Analog Output Range <sup>3</sup>	0.5 ~ 4.5	Vdc	
Response Time <sup>4</sup>	5	ms	
Pressure Sensor⁵			
Pressure Type	Gauge		
Pressure Range <sup>6</sup>	-5 ~ +40 / -5 ~ +100	cmH₂O	
Accuracy <sup>7</sup>	± 1.0	%FS	
Output <sup>2</sup>	Linear, Analog/l²C	ms	
Analog Output Range 8	0.5 ~ 4.5	Vdc	
Response Time <sup>4</sup>	5	ms	
Other Specifiactions			
Gas Type	Air		
Supply Voltage	5.0 ± 5%	Vdc	
Output Pin Shape 9	6 Pins, NS-TECH CD R-6		
Accuracy Compensated Temp.	-5 ~ <b>+</b> 65	°C	
Accuracy Compensated Altitude	-400 ~ +3000 (700 ~ 1060)	m (hPa)	
Storage Temp.	-40 ~ +85	°C	
Maximum Pressure 10	±30 / ±80	kPa	
Humidity	(0~100)%RH (No icing or condersation)		
Warming Up Time	<50	ms	

- 1, Other flow ranges are selectable, see detailed in section 5.
- 2, 14 bit for I<sup>2</sup>C.
- 3, Corresponding with flow rate of reverse full scale flow rate to forward full scale flow rate.
- 4, Update time for I<sup>2</sup>C is up to 1ms.
- 5, Pressure sensor is selectable.
- 6, Other pressure ranges are selectable.
- 7, Accuracy specifications apply over operating conditions. This specification represents the total combination of non-linearity, hysteresis, zero and span shift, repeatability & temperature effects.
- 8, Corresponding with flow rate of reverse full scale pressure to forward full scale pressure.
- 9, The 0.5 m output cable (SN6-50) is usually shipped together with the sensor.
- 10, Maximum pressure base on the working pressure of the pressure sensor.

Pressure range:  $-5 \sim +40 \, \text{cmH}_2\text{O}$ , proof pressure:  $\pm 30 \, \text{kPa}$ , burst pressure:  $\pm 40 \, \text{kPa}$ .


Pressure range: -5 ~ + 100 cmH<sub>2</sub>O, proof pressure: ±80 kPa, burst pressure: ±100 kPa.




### 3. Flow Characteristics

Flow Rate (SLPM)	Typical Voltage (Vdc)
-250	0.500
-200	0.900
-150	1.300
-100	1.700
-50	2.100
0	2.500
+50	2.900
+100	3.300
+150	3.700
+200	4.100
+250	4.500

Pressure (cmH <sub>2</sub> O)	Typical Voltage (Vdc)		
-5	0.500		
0	0.944		
+5	1.389		
+10	1.833		
+15	2.278		
+20	2.722		
+25	3.167		
+30	3.611		
+35	4.056		
+40	4.500		

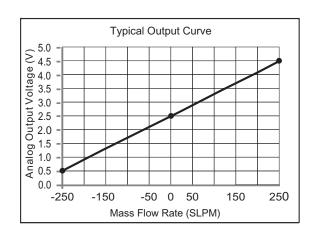

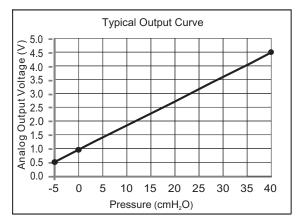

#### 4. Electrical Interface

#### 4.1 Pin Definition

The FS6122 provides a 6-pin interface. The output connecting cable comes with the sensor. The sensor pin layout is shown in Figure 4.1 and the cable color code is defined in Table 4.1.

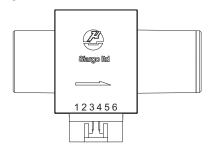


Figure 4.1: FS6122 pin layout.

## **4.2 Pin Description**

**VCC** and **GND**: The FS6122 requires a power supply of 5±5% Vdc. The voltage is internally filtered and regulated to power the circuit. The sensor consumes less than 15 mA normally and the minimum supply current must be larger than 15 mA.

**Pout**: The pressure analog output is selectable. The pin provides  $0.5 \sim 4.5$  Vdc corresponding with the specified pressure full scale range.

Table 4.1: FS6122 pin code.

	Pin	Color	Definition		
	1	White	Pout <sup>1</sup> , Pressure sensor analog outpu		
	2	Green	Fout, Flow sensor analog output		
	3	Black	GND, Ground		
_	4	Red	VCC, Power supply		
	5	Yellow	SCL (I <sup>2</sup> C)		
_	6	Blue	SDA (I <sup>2</sup> C)		

<sup>1,</sup> The pressure output is selectable.

**Fout**: The flow rate analog output. The pin provides 0.5 ~ 4.5 Vdc corresponding with the specified flow full scale range.


**SDA** and **SCL**: Serial data line and serial clock line.



## 5. Ordering Guide

#### **5.1 Sensor Selection**

The sensor part number is composed of the product model number and suffix indicating the full scale flow rate, mechanical connection, output format as well as the application gas. Refer the following for details.

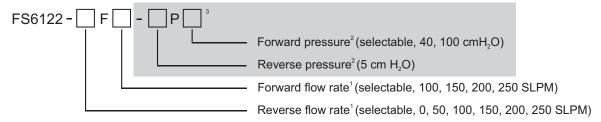

- 1 Max. flow rate number only, for example, 250 meaning full scale flow rate of 250SLPM. When reverse flow rate equal to 0, the sensor will be an uni-directional sensor.
- 2 Max. pressure number only, for example, 40 meaning full scale flow rate of 40 cmH2O. Other flow ranges are selectable
- 3 Pressure output is selectable, null for non-pressure models.
- 4 For instance, FS6122-50F200-5P40, means flow rate of -50 to +200 SLPM, pressure of -5 to +40 cmH2O.


#### 5.2 Order Contact and Customer Support

The sales offices are listed at the end of this document. For small quantities, the order can be placed either through Siargo website: www.siargo.com or the convenient sales office. For large quantities, please contact the sales office, distributors or sales representatives.

Siargo is making every effort to ensure the quality of the products. In case of questions and/or product supports, please contact customer service listed at the end of the document.

#### 6. Mechanical Dimensions


The FS6122 provides standard ISO-22mm medical connection and can be readily applied to CPAP machines. The sensor has a total size of 67 X 47 X 31 mm.

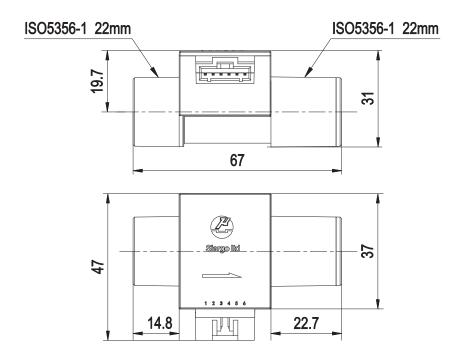




#### 7. I<sup>2</sup>C Communication

#### 7.1 I<sup>2</sup>C Connection

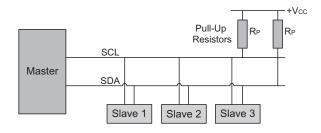

## 7.2 I<sup>2</sup>C Read and Write Sequences

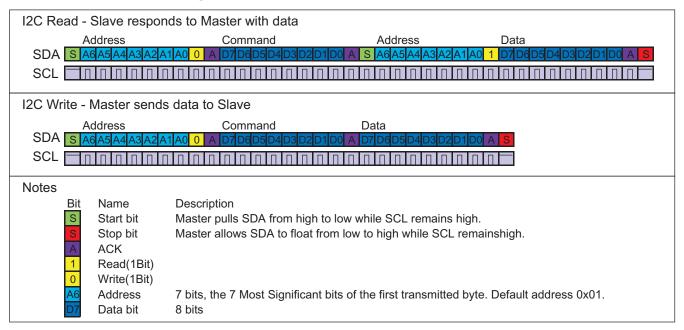

## 7.3 I<sup>2</sup>C Commands description

Command Byte (Hex)	Length	Command Name	Read <sup>1</sup>	Write <sup>2</sup>	Notes
02H	12	Read/Write the SensorSN		•	ASCII
03H	4	Readflow rate	0		Int32/1000 flow rate unit
23H	4	Readpressure	0		Int32/1000 pressure unit
04H	8	Readflow rate and pressure	<b>©</b>		Int32/1000 flow rate unit , and Int32/1000 pressure unit
05H	1	Read/Write the I <sup>2</sup> Caddress	0	0	Bit7 ~ Bit1 can be set <sup>3</sup>
0BH	1	Read/Write the filter depth	0	0	Int8, 0~254
1CH	1	Calibration the offset of flow rate		0	1byte, any data, NO flow in the pipe
24H	1	Calibration the offset of flow pressure		0	1byte, any data, NO flow in the pipe

- 1) O means the commands is available for any user.
- 2) means the commands is available only for the manufacturer.
- 3) The address is set with Bit7~Bit1, For instance, 08H means address of 4.
- 4) When reading information, add 80H to the command byte. For instance, Read flow rate command is 83H.



## **Important Notices**

#### **Wetted Materials and Compatibility**

The sensor body is made of medical compatible plastics. The sensor chip comprises of silicon, silicon nitride and silicon dioxide and the sensor chip surfaces are passivated with silicon nitride and silicon dioxide. The electronic sealing is provided by RTV (room temperature vulcanizing) silicone sealant WR-704 composed of HOCH<sub>3</sub> (SiO) CH<sub>3</sub>H.

## Cautions for Handling and Installations

The product at the time of shipment is fully inspected for product quality and meets all safety requirements. Additional safety measures during handling and installation should be applied. To prevent ESD (electrostatic discharge) damage and /or degradation, take customary and statutory ESD precautions when handling. Do power the product with the correct polarity, voltage & amperage. All precautions and measures for electrical voltage handling must apply. The product sealing is ensured to work under working pressure of 30 kPa and is leakage proof before the shipment. But cautions and further leakage test are important at installation as well since any leakage could cause severe safety issue.

This product contains no user serviceable components. Do not attempt to disassemble, substitute parts or perform unauthorized modifications to the product. Doing so will forfeit the terms of the warranty and cause the liability to any damages thereafter. It should only be serviced by authorized personnel. Upon requests, Siargo will provide necessary technical support and/or training of the personnel.

### **Cautions for Product Applications**

The product is designed for use with general purpose gases such as air and nitrogen. It is advised that the products are best used for non-explosive clean gases. The sensors cannot be used for gas metrology of fluoride or fluoride-containing gases. For updates of the product certification information, please contact the manufacturer. Use for other gases such as extreme corrosive and toxic may cause the product malfunctioning or even severe damages.

Don't expose the product's electronics other than the inner flow channel to any liquids, the unit does not have a water proof electronics. For medical sterilization procedure, please consult the manufacturer. Don't flow gas in conditions that can cause condensing water vapor to be trapped inside the unit during operation as the accuracy could be significantly influenced.

It is suggested to design your application so that nominal flow rate is approximately 70% of the full scale

flow rating of the sensor. Don't use a sensor with a flow range at the extreme cases, for instance, don't use a 200 SLPM sensor for a 2 SLPM application.

#### **Warranty and Liability**

(Effective January 2010)

Siargo warrants the products sold hereunder, properly used and properly installed under normal circumstances and service as described in this user manual, shall be free from faulty materials or workmanship for 180 days for OEM products, and 365 days for non-OEM products from the date of shipment. This warranty period is inclusive of any statutory warranty. Any repair or replacement serviced product shall bear the same terms in this warranty.

Siargo makes no other warranty, express or implied and assumes no liability for any special or incidental damage or charges, including but not limited to any damages or charges due to installation, dismantling, reinstallation or any other consequential or indirect damages of any kind. To the extent permitted by law. the exclusive remedy of the user or purchaser, and the limit of Siargo's liability for any and all losses, injuries or damages concerning the products including claims based on contract, negligence, tort, strictly liability or otherwise shall be the return of products to Siargo, and upon verification of Siargo to prove to be defective, at its sole option, to refund, repair or replacement of the products. No action, regardless of form, may be brought against Siargo more than 365 days after a cause of action has accrued. The products returned under warranty to Siargo shall be at user or purchaser's risk of loss, and will be returned, if at all, at Siargo's risk of loss. Purchasers or users are deemed to have accepted this limitation of warranty and liability. which contains the complete and exclusive limited warranty of Siargo, and it shall not be amended, modified or its terms waived except by Siargo's sole

This warranty is subject to the following exclusions:

- (1) Products that have been altered, modified or have been subject to unusual physical or electrical circumstances indicated but not limited to those stated in this document or any other actions which cannot be deemed as proper use of the products;
- (2) Siargo does not provide any warranty on finished goods manufactured by others. Only the original manufacturer's warranty applies;
- (3) Products re-sold to the third parties.



#### **Contact Information**

#### **Headquarters**

Siargo Ltd. 3100 De La Cruz Boulevard, Suite 210, Santa Clara, California 95054, USA Phone: +1-408-969-0368

Email: info@Siargo.com

#### Representative in US

Servoflo Corporation 75 Allen Street, Lexington, MA 02421, USA Phone: +1-781-862-9572 Email: info@servoflo.com

#### Representative in Japan

Marubeni Information Systems Co., Ltd.
Device Solutions Department
14th Floor, Shinjuku Garden Tower,
3-8-2, Okubo, Shinjuku-ku, Tokyo 169-0072, Japan

Phone: +81-3-4243-4160 Fax: +81-3-4243-4198

Email: siargo@marubeni-sys.com

#### **Representative in Europe**

IDENTIC Stefanie Trojan In der Siedlerruh 24 69123 Heidelberg / Germany Phone: +49-(0)6221-7509777 Fax: +49-(0)6221-7509779

Email: info@identic.de

This document is intended for distribution to customer(s) only. Siargo and its subsidiaries reserve the rights to change the specifications and/or descriptions without prior notice.